



Rotaid **PLUS** Indoor Cabinet for Defibtech



English
version

Product sheet

The Rotaaid PLUS Indoor Cabinet for Defibtech provides a robust, durable, and recognizable indoor storage solution for all of your Lifeline AEDs. The characteristic round shape makes your life-saving equipment easy to find, and just a quick twist is needed to gain full access to the AED. The green tinted polycarbonate cover comes with the international AED logo and the directive arrows and signage provide user guidance. The cabinet comes with a battery powered, audible alarm system to alert bystanders when the system is in use.

The unique cover is equipped with recesses for optimal grip with distinctive yellow striping and Defibtech's logo. Its solid construction and UV filters ensure a long-lasting appearance and provide sturdy protection for your life-saving equipment. The translucent material makes sure your Lifeline AED and its Active Status Indicator (ASI) are visible at all times.



The cabinet is opened by a quick and intuitive twist that removes the cover from the backplate. A nylon strap holds the cover to the backplate, and the responder has full access to the AED. The AED inside is safely positioned on a metal stand and can be grabbed in an instant. All models of the Lifeline AED series, and the Lifeline VIEW series can be used or stored with or without their carrying cases. A solid steel hook, attached to the bottom of the cabinet, can hold a tag, seal, or other rescue accessory.

Modern production techniques, and the use of recycled materials, allow the cabinet to meet the highest requirements when it comes to sustainability. Because of the battery powered alarm system, a power socket is not required, enabling an easy installation process.



Specifications

Application: Indoor use

WARNING: not for outdoor use

Dimensions

- Length: 47 cm 18,5"
- Width: 47 cm 18,5"
- Depth: 18 cm 7,1"
- Weight: 1,7 kg 3,75 lb

Backplate

- 100% recycled ABS (acrylonitrile butadiene styrene), 4 mm thick
- Solid steel hook on the bottom to install tamper seal, instructions or other accessories
- A metal stand allows for safe placement of your AED

Alarm system

- The pre-installed audible alarm system, triggered by a coin sized, low power magnet in the cover, has an ABS housing and is affixed to the backplate
- The unit is powered by two AA-batteries (included)
- On the bottom of the alarm housing, a red button can be found to temporarily pause the alarm if required
- The alarm loudness level can be selected via a white switch on the bottom of the alarm housing
Set to the left, produces an 84 dB level
Set to the right produces a 98 dB level

International warranty

- 12 months

Packaging

- Box size 50 x 50 x 25 centimeters / 19,7 x 19,7 x 9,8 inches
- All fastening materials included
- One cabinet per box

Cover


- 100% recycled green tinted polycarbonate, 5 mm thick, fully transparent, fitted with UV filters for long-lasting transparency and UV blocking
- A nylon strap is attached to the cover and connected to the backplate, and is fitted with brass rings
- 4 direction arrows and sticker provide guidance for opening the device
- 6 indents provide grip for fast and easy opening

Compatible Defibtech AEDs

- Lifeline AED / Lifeline AUTO
- Lifeline VIEW / Lifeline VIEW AUTO
- Lifeline ECG
- Lifeline PRO

Compatible Defibtech Carrying Cases:

- Soft Carrying Case for DDU-100 series
- Soft Carrying Case for DDU-2000 series
- Semi-Rigid Carrying Case for DDU-2000 series

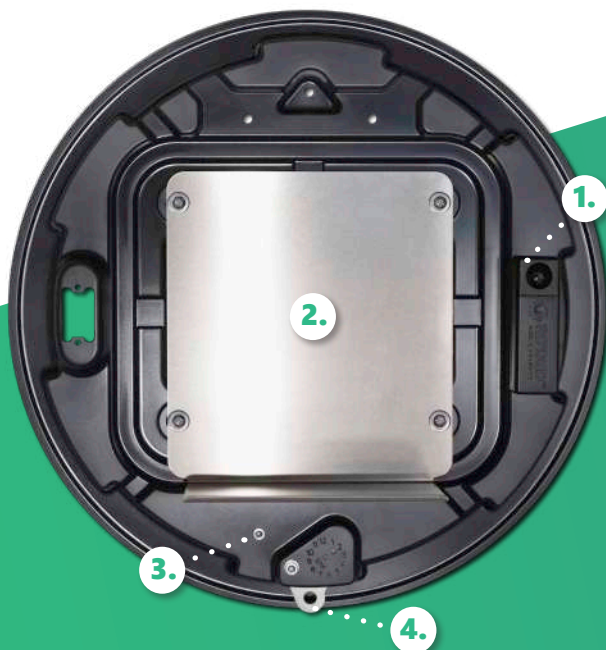

1. Alarm system
2. Metal stand
3. Cover strap holder
4. Stainless steel hook

Placement of Defibtech Lifeline AEDs







Lifeline AED



Lifeline VIEW





 rotaid.com  info@rotaid.com  +31 (0)478 769 052
 Keizersveld 52D 5803 AN, Venray The Netherlands

THIS DOCUMENT IS FOR INTERNATIONAL DISTRIBUTION ONLY

MKT-GD052-100-EN
ELECTRONIC DISTRIBUTION